

Fiber Bundle Cleaning




- 1. Click Menu Icon.
- 2. Click Scanner Control (Password: fireadmin).

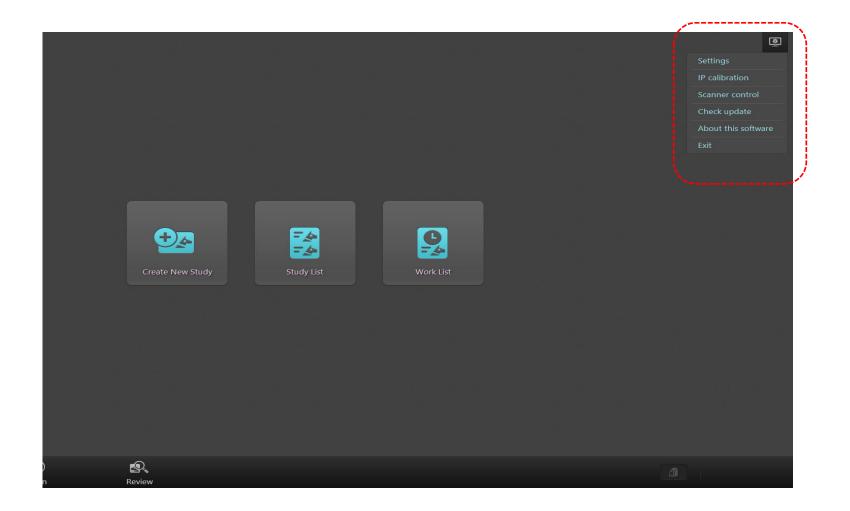


- 3. Click Forward.
- 4. Click Stop after 10 [sec].

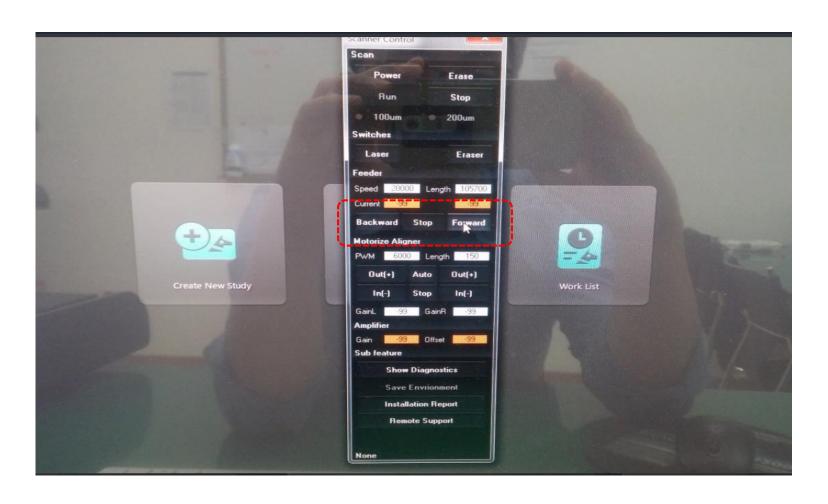

- 5. Turn off power
- 6. Remove the Top Cover and Internal Cover.
- 7. Disconnect three cassette lock cables from the cassette lock board.
- 8. Remove four screws from Plexi Guide.
- 9. Carefully lift up the Plexi Guide and then remove the RFID cable on the RFID board.

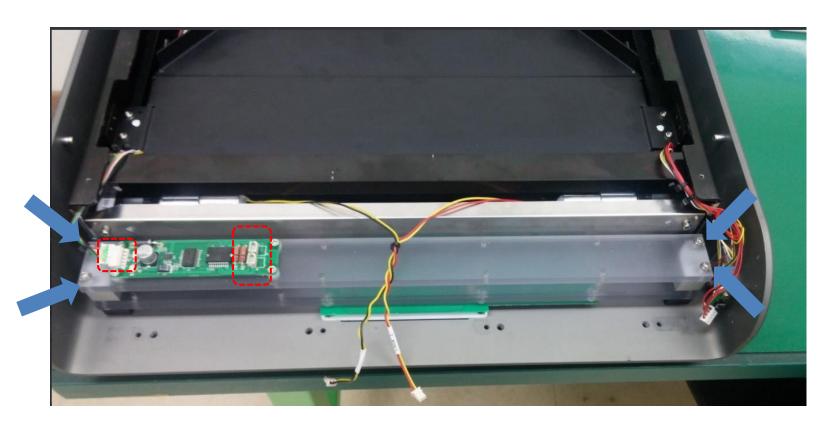

- 10. Remove two bolts.
- 11. Remove shield plate.

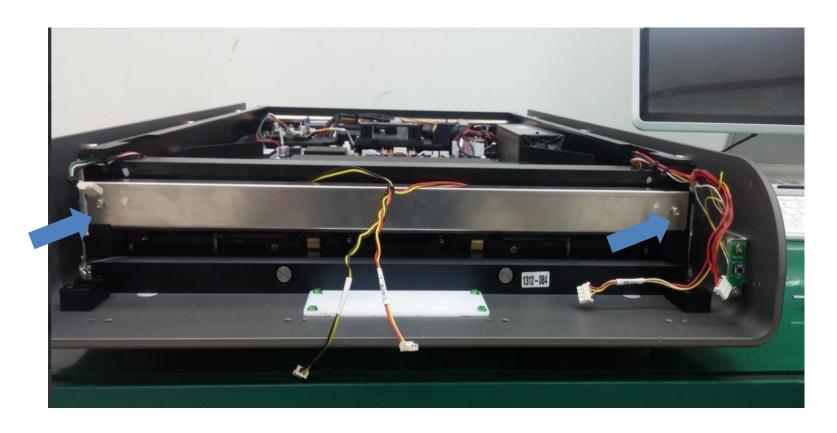


- 12. Clean the fiber bundle.
- 13. Clean it around fiber bundle.

Note) Should use the soft tissue with Alcohol when you clean the fiber bundle.



14. Clean the fiber bundle cover.

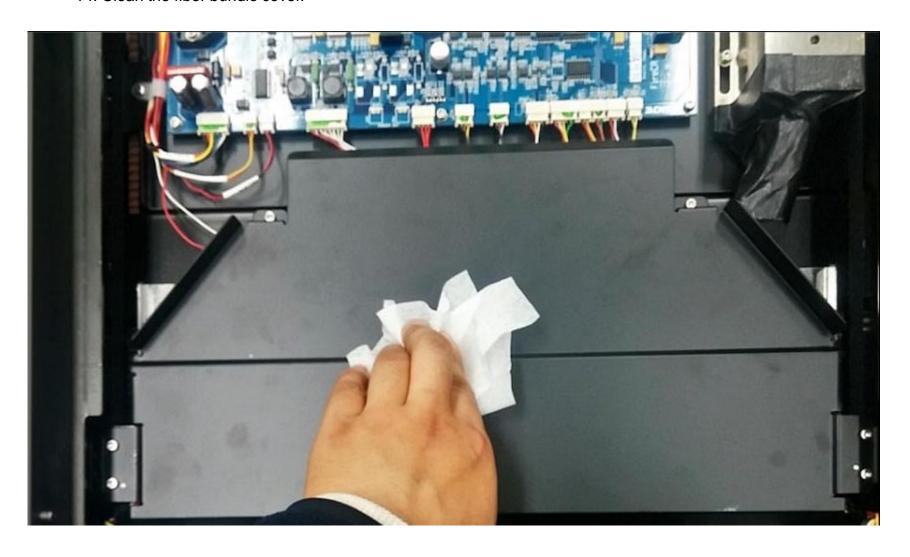

- 15. Clean the reflection mirror.
- 16. Clean F-theta lens.

