

QUANTORDENT IMAGING SOFTWARE

- ☑ Fast, high-definition images
- ☑ Easy image viewing and manipulation
- ☑ Free software upgrades



Easy-to-Use Imaging Software

With the QuantorDent Imaging Software image acquisition, processing, and management have never been easier. The easy-to-use software, network/TWAIN interface and free software upgrades makes it a safe investment. The QuantorDent Imaging Software is the software for 3DISC's FireCR Dental PSP Readers.

Easy Image Acquisition & Processing

Experience the high diagnostic value of your FireCR Dental PSP Reader. After just a few seconds of scanning, the image appears. The chronological and intuitive workflow helps you crop, enhance, increase brightness and contrast, and perform other adjustments on the images to enable a more precise diagnosis. Regions of interest can be highlighted, annotations added and images can easily be rearranged or deleted.

Intuitive Interface & Study Management

The intuitive user interface allows simple integration into your clinic's existing patient management system - providing the full range of capabilities needed to improve your productivity and workflow. The QuantorDent Imaging Software allows clinics to quickly input patient data—or access it directly from their patient management system.

Download patient and exam information from a work list, which can be automatically matched, eliminating the need for data entry. Manually input patient and exam information, if needed.

Embedded Viewer for Immediate Diagnosis

The embedded viewer is your everyday viewing tool, allowing for immediate diagnosis. The viewer contains tools such as annotations, contrast, brightness, measurements, drawings. For smaller facilities without a HIS, RIS, or PACS, a local database enables direct management of patients and studies. QuantorDent Imaging Software is available in a variety of languages.

☑ HIGHLIGHTS

- Fast, high-quality image processing designed for CR
- Streamlined intuitive workflow
- Full image manipulation functionality
- Add and customize dental templates
- Fast and easy image verification
- Simple data management
- Interfaces with existing patient management system
- Secure, high-speed local database for facilities without HIS, RIS or PACS
- Full range of image distribution & output options
- Touch-screen user interface
- Automation using DICOM work list server
- Support various screen resolutions for notebook or tablet PCs

Quantor Database Options

Use your own standalone database or save your studies to a shared database with QuantorServer and view them with QuantorDent Lite.


The QuantorDent Lite software is installed on any PC wanting to access the QuantorServer shared database on the network. A license dongle is not required.

Once you run QuantorDent Lite it will automatically search the network for an available QuantorServer and set itself up.

A Standalone Database

The local database can only be accessed from the PC on which the database is stored. QuantorDent must be installed on each PC.


- No access to databases or sharing of images across networks
- Each database is separately stored locally on PC
- *LDB = Local Database

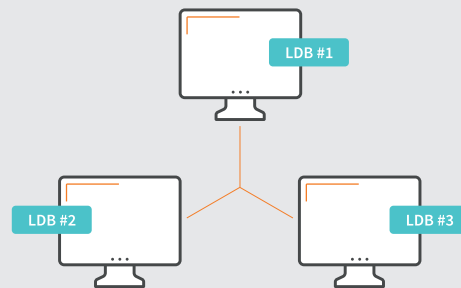

B Shared Database

Save your studies to a shared database with QuantorServer and view them with QuantorDent Lite. Allows you to share studies between all workstations on a network.

- All PCs can access the same database
- Operating PCs run QuantorDent
- Viewing-only PCs run QuantorDent Lite

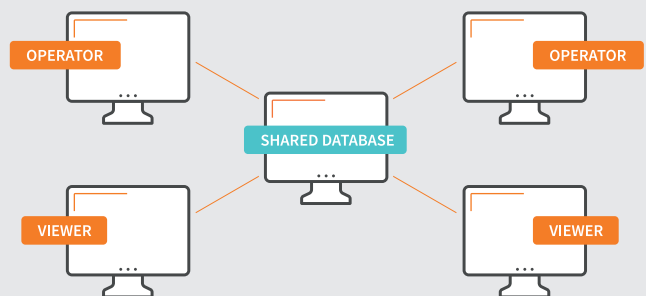

FireCR, FireCR+, FireCR Flash, FireCR Spark, FireCR Dental, FireID, Quantor, QuantorMed, QuantorMed+, QuantorVet, QuantorVet+, QuantorDent, QuantorDentVet, QuantorDB & QuantorShare are trademarks of Digiray Corp.

We reserve the right to modify specifications without prior notice. © 2019 Copyright Digiray.